

SepraPor®

Cartridge Insertion Tool


Cartridge Insertion Tools are available to facilitate insertion of either a 4" or 6" diameter Seprapor® hollow fiber cartridge into a stainless steel housing. These tools ensure an easy, ergonomic installation of the filter cartridge into the housing.

Features & Benefits

- Optional tools simplify the installation of either a 4" or 6" cartridge into a stainless steel housing
- Applies even pressure across cartridge to prevent crooked installation and o-ring displacement during insertion
- Viewing holes allow visual monitoring for proper insertion and application of lubricating solution (e.g. 70% isopropanol)
- Eliminates operator fatigue from manually pushing cartridge into place
- Durable sanitary 316L stainless steel design resists corrosion and allows for easy cleaning

Materials of Construction

Body:	Stainless Steel
Threaded Rod:	Stainless Steel
Center nut:	Brass
Handle:	Stainless Steel
Handle Grips:	Rubber

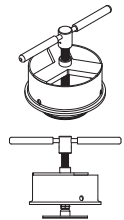

Ordering Information

HX4001-00	4" Cartridge Seating Tool, 316 Stainless Steel
HX6001-01	6" Cartridge Seating Tool, 316 Stainless Steel

